

Articulating Your Way!



MainStay™ Disposable Articulator System


Whip Mix Corporation has raised the bar on disposable articulation with the Mainstay™ Articulator System. This new line of articulators was designed to save time, money and waste; and does that with the incomparable stability of the MainStay™ Pins.

The Mainstay™ Disposable Articulator System does not require model formers or a gypsum base, and adjusts to accommodate any occlusal relationship. Its base was designed without external walls or raised edges to prevent interference with the seating of dies and model segments. This unique feature also means that it can easily be used with gypsums of all expansion ratios.

Choose between Pin Models in Full Arch, Anterior and Quadrant sizes. Pinless Model offered in a Quadrant size only.


MainStay™ Disposable Articulator System

The benefits are easy to Articulate!

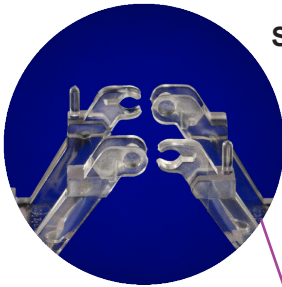


Strong Hinges

- Easy to separate and assemble
- Durable- resists breakage
- Can be heated to reposition occlusal relationship

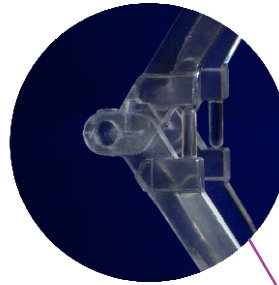


Built-in Support Rods

- Provide a vertical stop, ideal for edentulous cases

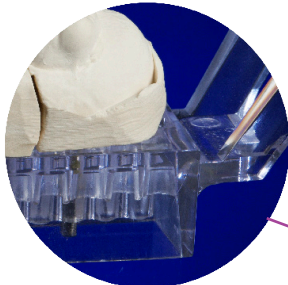


No External Walls

- Flat tray for easy removal and replacement of segments
- A feature in both pin and pinless designs

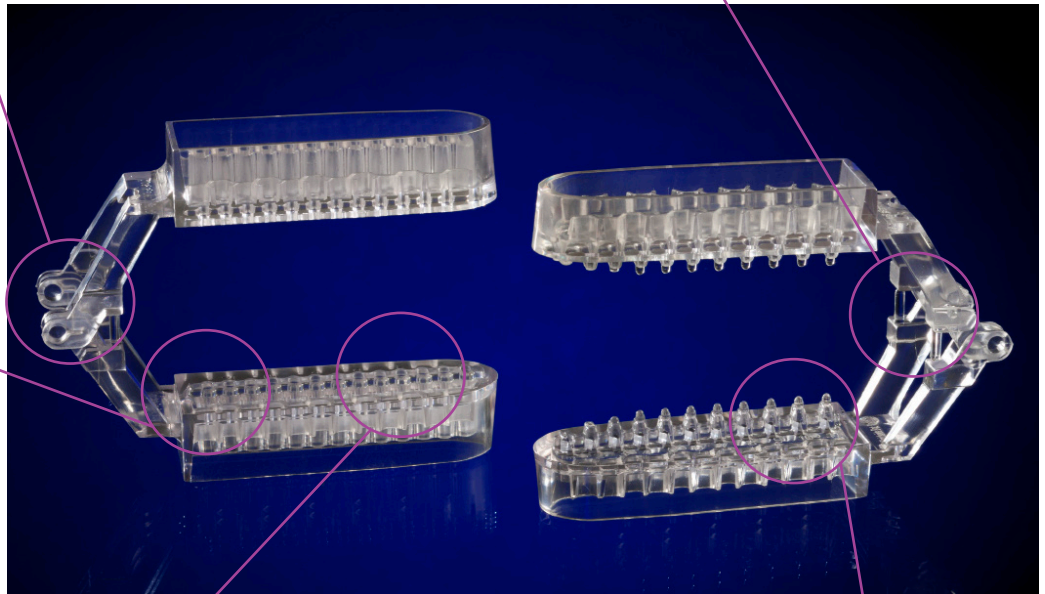


MainStay™ Pins

Unique tapered blade shape maintains precise position of dies

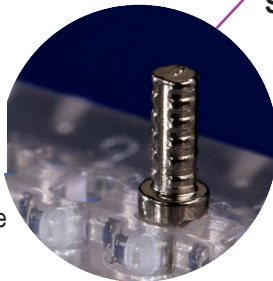

Stable Fit with Pin or Pinless System

- Raised anti-rotational spine prevents movement of dies
- Numbered bases for quick reference of die position
- Excellent for scanning
- Pinless design eliminates pinning, drilling, gluing
- Combine Pin Quadrant for working model and Pinless Quadrant for opposing

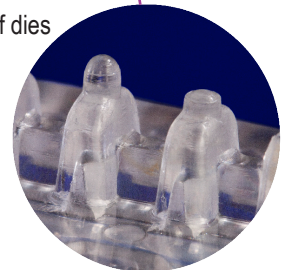

Full Arch Pin Item # 08251

Set Includes: 25 Articulators, 100 MainStay™ Pins and 3 Support Rods

Anterior Pin Item # 08253

Set Includes: 25 Articulators, 100 MainStay™ Pins and 3 Support Rods

Quadrant Pin Item # 08252

Set Includes: 50 Articulators, 100 MainStay™ Pins and 3 Support Rods

Quadrant Pinless Item # 08254

Set Includes: 50 Articulators and 3 Support Rods

Sold Separately:



Support Rods

Item # 08255, set of 20
Length: 1-7/8"

MainStay™ Pins

Item # 08256, set of 1000

To receive a Free Sample Box of MainStay™ Disposable Articulators, call **800-626-5651 x1270** or visit **whipmix.com**



more reasons. one source

Whip Mix Corporation • 361 Farmington Ave • Louisville, KY 40209 • 502-637-1451 • 800-626-5651 • www.whipmix.com